

-10°C~-25°C Biomedical Freezer

New Upgrade

Stable Performance • High Quality



DW - YL90



Brand New Appearance, the Second Generation of High-quality Product

The polyurethane foam door has great thermal insulation performance.




Multiple alarms ensure safety of storage

The freezer has multiple audible & visual alarm functions, such as high temperature/low temperature alarm, door opening/closing alarm, sensor failure alarm, etc. and can guarantee security of samples in a comprehensive way.



Temperature Control

The high-precision computerized temperature control system ensures an adjustable temperature within a range from -10 to -25°C inside the cabinet.

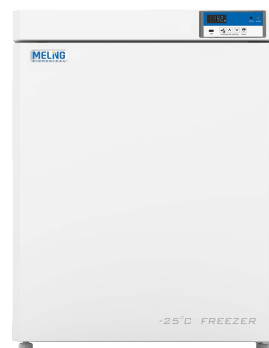



Control Panel

The high-definition digital temperature display can realize a temperature display precision of 0.1°C and is convenient for observation;



Alarm Upgrading

The high temperature/low temperature alarm, door opening/closing alarm, sensor failure alarm and power outage alarm can ensure security of samples comprehensively;



Refrigeration System

The compressor supplied by a famous brand and high-efficiency refrigerant R600a can reduce energy consumption to a large extent;




Human-oriented

The invisible door handle can prevent damage during transportation and save space;
The 3 drawers are convenient for separate storage of different types of things and can prevent cross influence among various things.


✓ Scope of Application

Suitable for freezing of ice bars and storage of various things needing refrigerated storage such as blood plasma, reagent, etc. Suitable for use in hospitals, clinics, health and disease prevention systems, blood banks, laboratories in colleges & universities, the frozen food industry and catering industry, etc.

🔗 Performance Data / Cooling Curve

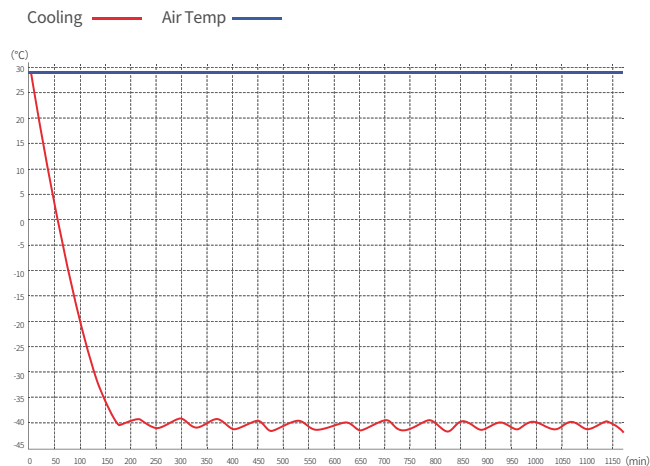

📐 External Dimensions

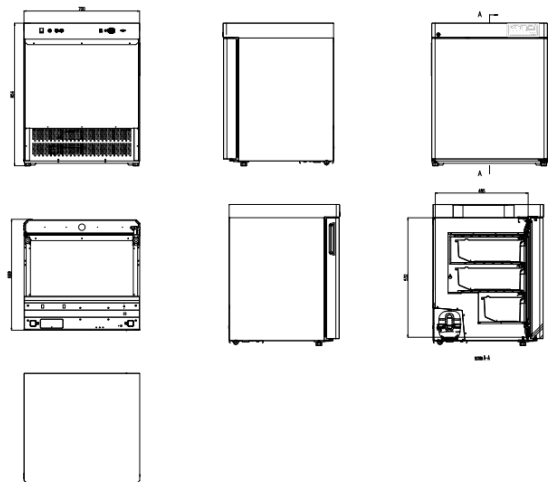

-25°C Low-temperature Medical Freezer	
Model	DW-YL90
Cabinet Type	Upright
Capacity(L)	90
Internal Size(W*D*H)mm	500*485*532
External Size(W*D*H)mm	700*669*854
Package Size(W*D*H)mm	765*735*935
NW/GW(Kgs)	54/67
Performance	
Temperature Range	-10~-25°C
Ambient Temperature	16-32°C
Cooling Performance	-25°C
Climate Class	N
Controller	Microprocessor
Display	Direct display
Refrigeration	
Compressor	1
Cooling Method	Direct Cooling
Defrost Mode	Manual
Refrigerant	R290
Insulation Thickness(mm)	100
Construction	
External Material	PCM
Inner Material	Aluminum plate with spraying
Drawers	3(ABS)
Door Lock with Key	Yes
Access Port	1
Castors	2+(2 leveling feet)
Alarm	
Temperature	High/low temperature,High ambient temperature
System	Power failure, Low battery
Construction	Sensor failure,Door ajar,Built-in USB datalogger failure,Communication failure
Electrical	
Power Supply(V/HZ)	220-240~/50
Rated Power(W)	80

*The model, parameters and performance specified in this brochure may be changed without prior notice because of product upgrading.

*There may be differences between the product images shown in this brochure and the actual products. When you are buying any product, please check the actual product.